

Dr. Fritz®

**Scissors, dissecting
and grasping forceps**
Minimal invasive surgery

Inst 1-03E



- Standard handles
- easy dismantable
- long life span - one system
- reusable

The modular system of grasping forceps and scissors provides the user with a range of products which are of lasting value and that satisfy every requirement as regards handling, cleaning and preparation.

The system

The modular construction makes it possible to assemble an optimal instrument for surgery, giving a choice of different jaw inserts, tubular sheaths and handles.


Thanks to the accurately manufactured bayonet lock between the jaw insert and the sheath, as well as "twist connections" between the sheath and the various handles, the instruments can be dismantled quickly. This makes it possible to clean parts of the equipment that would otherwise be difficult to access, and also makes procurement of individual modules cost-effective.

Sheaths and jaw inserts of the instruments can be rotated through 360°.

The uncomplicated construction and the use of state of the art materials guarantee a long service life of the instruments.

Jaw inserts

- 5 and 10 mm diameter
- working lengths of 23 cm, 31 cm and 43 cm

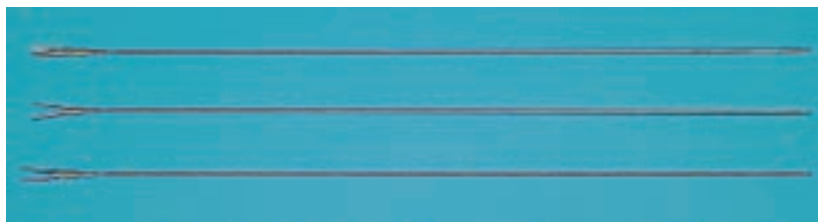

Tubular sheaths

- 5 mm und 10 mm diameter
- working lengths of 23 cm, 31 cm and 43 cm
- rotateable through 360°, or fixed
- insulated or available in stainless steel



Handles

- without, with external or with internal ratchet
- with HF-connector
- colour coding for scissors



Scissors and biopsy forceps, 5 mm



Working length

31 cm

23 cm

43 cm

METZENBAUM scissors

curved "Super Cut", HF-connector

Complete instrument

Insert
Sheath
Handle

531.01 H 05	521.01 H 05	541.01 H 05
53 M 01	52 M 01	54 M 01
53 R 1	52 R 1	54 R 1
H 05	H 05	H 05

MINI-METZENBAUM scissors

curved "Super Cut", HF-connector

Complete

Insert
Sheath
Handle

531.02 H 05	521.02 H 05	541.02 H 05
53 M 02	52 M 02	54 M 02
53 R 1	52 R 1	54 R 1
H 05	H 05	H 05

Hook scissors

straight, HF-connector

Complete

Insert
Sheath
Handle

531.04 H 05	521.04 H 05	541.04 H 05
53 M 04	52 M 04	54 M 04
53 R 1	52 R 1	54 R 1
H 05	H 05	H 05

Scissors

curved, Tungsten Carbide

Complete

Insert
Sheath
Handle

531.03 H 05	521.03 H 05	541.03 H 05
53 M 03	52 M 03	54 M 03
53 R 1	52 R 1	54 R 1
H 05	H 05	H 05

Scissors

straight, Tungsten Carbide

Complete

Insert
Sheath
Handle

531.06 H 05	521.06 H 05	541.06 H 05
53 M 06	52 M 06	54 M 06
53 R 1	52 R 1	54 R 1
H 05	H 05	H 05

MINI scissors, 3/5 mm

curved

Complete

Insert
Sheath
Handle

531.08 H 05	521.08 H 05	541.08 H 05
53 M 08	52 M 08	54 M 08
53 R 1	52 R 1	54 R 1
H 05	H 05	H 05

Biopsy forceps

without spikes

Complete

Insert
Sheath
Handle

531.60 H 01	521.60 H 01	541.60 H 01
53 M 60	52 M 60	54 M 60
53 R 1	52 R 1	54 R 1
H 01	H 01	H 01

Biopsy forceps

with spikes

Complete

Insert
Sheath
Handle

531.61 H 01	521.61 H 01	541.61 H 01
53 M 61	52 M 61	54 M 61
53 R 1	52 R 1	54 R 1
H 01	H 01	H 01

Please note!

Due to the different jaw types the working length of the complete instrument can be different to other types

Dissecting forceps and grasping forceps, 5 mm



Working length

31 cm

23 cm

43 cm

Dissecting forceps

short, curved

Complete instrument

Insert
Sheath
Handle

531.40 H 01	521.40 H 01	541.40 H 01
53 M 40	52 M 40	54 M 40
53 R 1	52 R 1	54 R 1
H 01	H 01	H 01



Dissecting forceps, PETELIN

short, curved, only to use with trocar cannulas with an inner diameter of 8 mm

Complete

Insert
Sheath
Handle

531.37H 01	521.37 H 01	541.37 H 01
53 M 37	52 M 37	54 M 37
53 R 1	52 R 1	54 R 1
H 01	H 01	H 01

Dissecting forceps

straight

Complete

Insert
Sheath
Handle

531.42 H 01	521.42 H 01	541.42 H 01
53 M 42	52 M 42	54 M 42
53 R 1	52 R 1	54 R 1
H 01	H 01	H 01

Grasping forceps

small front hooks, with ratchet

Complete

Insert
Sheath
Handle

531.20 H 21	521.20 H 21	541.20 H 21
53 M 20	52 M 20	54 M 20
53 R 1	52 R 1	54 R 1
H 21	H 21	H 21

Grasping forceps

two pairs of hooks, with ratchet

Complete

Insert
Sheath
Handle

531.21 H 21	521.21 H 21	541.21 H 21
53 M 21	52 M 21	54 M 21
53 R 1	52 R 1	54 R 1
H 21	H 21	H 21

Grasping forceps

universal, with ratchet

Complete

Insert
Sheath
Handle

531.22 H 11	521.22 H 21	541.22 H 21
53 M 22	52 M 22	54 M 22
53 R 1	52 R 1	54 R 1
H 11	H 21	H 21

Grasping forceps

straight, tapered, with ratchet

Complete

Insert
Sheath
Handle

531.24 H 21	521.24 H 21	541.24 H 21
53 M 24	52 M 24	54 M 24
53 R 1	52 R 1	54 R 1
H 21	H 21	H 21

Claw grasping forceps

with ratchet

Complete

Insert
Sheath
Handle

531.30 H 21	521.30 H 21	541.30 H 21
53 M 30	52 M 30	54 M 30
53 R 1	52 R 1	54 R 1
H 21	H 21	H 21

Dissecting forceps and grasping forceps, 5 mm



Working length

31 cm

23 cm

43 cm

Complete instrument

Insert
Sheath
Handle

Atraumatic grasping forceps

curved, with ratchet

531.23 H 21	521.23 H 21	541.23 H 21
53 M 23	52 M 23	54 M 23
53 R 1	52 R 1	54 R 1
H 21	H 21	H 21

Complete

Insert
Sheath
Handle

Atraumatic grasping forceps

straight, with ratchet

531.25 H 21	521.25 H 21	541.25 H 21
53 M 25	52 M 25	54 M 25
53 R 1	52 R 1	54 R 1
H 21	H 21	H 21

Complete

Insert
Sheath
Handle

Grasping forceps (BABCOCK)

with ratchet

531.27 H 22	521.27 H 22	541.27 H 22
53 M 27	52 M 27	54 M 27
53 R 1	52 R 1	54 R 1
H 22	H 22	H 22

Complete

Insert
Sheath
Handle

Grasping forceps (DEBAKEY)

with ratchet

531.35 H 21	521.35 H 21	541.35 H 21
53 M 35	52 M 35	54 M 35
53 R 1	52 R 1	54 R 1
H 21	H 21	H 21

Complete

Insert
Sheath
Handle

Bowel grasping forceps

without HF-connector, with ratchet

531.36 H 22	521.36 H 22	541.36 H 22
53 M 36	52 M 36	54 M 36
53 R 1	52 R 1	54 R 1
H 22	H 22	H 22

Complete

Insert
Sheath
Handle

Oviduct grasping forceps

with ratchet

531.44 H 21	521.44 H 21	541.44 H 21
53 M 44	52 M 44	54 M 44
53 R 1	52 R 1	54 R 1
H 21	H 21	H 21

Complete

Insert
Sheath
Handle

Mesometrium grasping forceps

HF-connector, with ratchet

531.45 H 21	521.45 H 21	541.45 H 21
53 M 45	52 M 45	54 M 45
53 R 1	52 R 1	54 R 1
H 21	H 21	H 21

Complete


Insert
Sheath
Handle

Adhesion and dissecting forceps

with ratchet

531.46 H 21	521.46 H 21	541.46 H 21
53 M 46	52 M 46	54 M 46
53 R 1	52 R 1	54 R 1
H 21	H 21	H 21

Scissors, dissecting forceps and grasping forceps, 10 mm

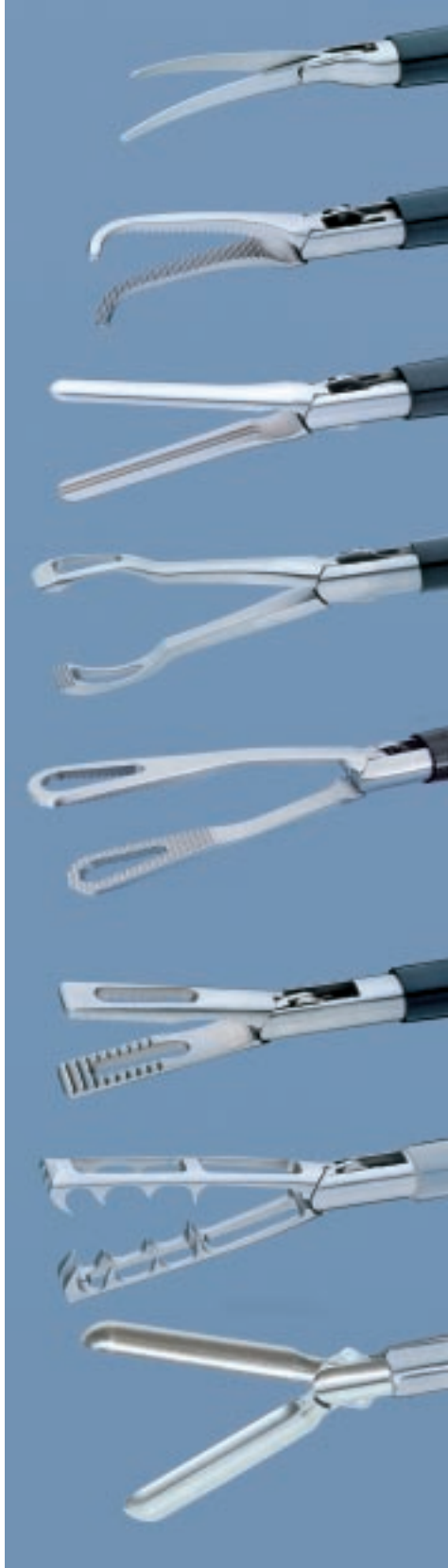

Working length

(total length incl. jaw)

27cm (31cm)

(standard version)

31 cm (36cm)

(long version)

MAYO scissor

curved, "Super Cut"

Complete instrument

Insert
Sheath
Handle

121.11 H 05

12 M 11

12 R 1

H 05

131.11 H 05

13 M 11

13 R 1

H 05

Dissecting forcep (OVERHOLD/MIXTER)

Complete

Insert
Sheath
Handle

121.48 H 01

12 M 48

12 R 1

H 01

131.48 H 01

13 M 48

13 R 1

H 01

Grasping forceps (DEBAKEY)

with ratchet

Complete

Insert
Sheath
Handle

121.31 H 21

12 M 31

12 R 1

H 21

131.31 H 21

13 M 31

13 R 1

H 21

Grasping forceps (BARCOCK)

with ratchet

Complete

Insert
Sheath
Handle

121.26 H 21

12 M 26

12 R 1

H 21

131.26 H 21

13 M 26

13 R 1

H 21

Atraumatic Grasping forceps

with ratchet

Complete

Insert
Sheath
Handle

121.49 H 21

12 M 49

12 R 1

H 21

131.49 H 21

13 M 49

13 R 1

H 21

Grasping forceps

universal, with ratchet

Complete

Insert
Sheath
Handle

121.25 H 21

12 M 25

12 R 1

H 21

131.25 H 21

13 M 25

13 R 1

H 21

Grasping forceps

non-rotating metal sheath, with ratchet

Complete

Insert
Sheath
Handle

121.32 H 21

12 M 32

12 R 1

H 21

131.32 H 21

13 M 32

13 R 1

H 21

Spoon forceps

large, with ratchet

Complete

Insert
Sheath
Handle

126.70 H 21

12 M 70

12 R 1

H 21

136.70 H 21

13 M 70

13 R 1

H 21



DE-3-456

Cleaning brush for long tubes, size 5-10 mm
pg. with 6 pcs.



Dr. Fritz GmbH · www.dr-fritz.com · info@dr-fritz.com
78532 Tuttlingen / Germany · Tel. +49/74 62/94 50 41 · Fax +49/74 62/94 50 45



Specifications and design are subjects to change without prior notice

07.10